PROBERG

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ENGLISH

Contents

Introduction	15
Important safety precautions	15
What's in the box	17
Parts of the product.....	17
Charging the Proberg MG-999.....	19
Instructions for use for the Proberg MG-999	19
Technical specifications for the Proberg MG-999.....	20
Cleaning and maintenance	20
Troubleshooting.....	21
Warranty policy.....	21
Electric and electronic devices	22
Guidance for safe removal of batteries & accumulators.....	23
Battery safety information.....	23
Customer service	24

Introduction

Dear customer,

Thank you for choosing Donnerberg Proberg MG-999 percussion massager with rechargeable battery. Please read the manual instructions carefully and save it for later use, be sure to make them accessible to other users. Please follow the instructions and use the product as described.

Proberg MG-999 is a battery-operated handheld massage device that offers deep tissue massage, which can help relieve pain and stiffness, reduce fatigue, enhance blood circulation, speed up muscle recovery, increase range of motion, improve flexibility. It can be used before and after a workout or physical activity. The device is only intended for domestic / private use, not for commercial and medical use. If you are unsure whether the massage gun is for you, please consult your doctor.

Important safety precautions

Before using electrical devices, basic safety precautions should always be followed, including the following:

READ ALL THE INSTRUCTIONS CAREFULLY BEFORE USE.

- The device can be used by adults only.
- This device should not be used by children, invalids or disabled persons without adult supervision.
- Children are not allowed to play with the device. Cleaning and user maintenance **MUST NOT** be done by children without adult supervision.
- After each use, the device must always be switched off and unplugged. Always unplug the device when not in use and before putting on and taking off the attachments.
- **DO NOT** reach for the massager that has fallen into water or other liquids. Keep dry – **DO NOT** operate in wet or moist conditions. **DO NOT** use while bathing or having a shower.
- **DO NOT** operate where aerosol (spray) products are being used or where oxygen is being administered.
- **DO NOT** use the massager while driving.
- **DO NOT** use the product before and/or while in bed. Massage has a stimulating effect and may delay sleep.
- **NEVER** cover the device while operating.
- The voltage and electrical to which the massager is connected must always comply with the laws in force.
- Use the Proberg MG-999 only on a dry surface of the body by lightly pressing and moving across targeted body areas. Use the Proberg MG-999 only on soft tissue making sure not to produce any discomfort or pain. **DO NOT** use on your head or bony part of the body.
- **DO NOT** use the massager on irritated, inflamed, swollen areas or skin eruptions. Stop if massage feels unpleasant or painful.

- NEVER use this device on or near your eyes and other sensitive areas.
- Keep long hair, clothing and jewellery away from the device while in use.
- Never leave the Proberg MG-999 operating or charging unattended, especially if children are present.
- DO NOT use any attachments not recommended or specifically provided by Donnerberg.
- The device has a heated surface. Persons insensitive to heat must be careful when using the appliance.
- NEVER use this massager if the battery or cable is damaged. If the massager is damaged, dropped, or broken, discontinue the further use.
- Do not attempt to repair the Proberg yourself, as any malfunction must be referred to the Donnerberg Service Centre.
- This device contains a non user-replaceable battery.
- This device contains a battery that is only replaceable by skilled persons or authorized service centres.
- For the purposes of recharging the battery, only use the power adapter provided with this device.
- Keep the device away from fire and excessive temperatures. DO NOT use the device in places where the temperature is above 40°.
- Extensive use may lead to the product's overheating and shorter life. If this occurs, please discontinue using the products and leave it to cool down before operating.
- DO NOT use the Proberg MG-999 continuously on one area if the massage feels unpleasant or uncomfortable, as this massager is designed to be used in a continuous motion. If bruising occurs, stop using the device.
- Use this device for the intended use only. DO NOT use this device to supplement or replace medical treatment.

You are strongly advised to consult your doctor before using the Proberg MG-999 if you:


- are pregnant;
- have a pacemaker, electronic implants, artificial joint(s) or limb(s);
- suffer from one or more of the following illnesses and complaints: circulation disorders, varicose veins, acute thrombosis, heart condition or cardiovascular disease, herniated disks, osteoporosis, spondylosis, epilepsy, severe migraines, diabetes or other sensory impairment, cancer;
- have had recent operation or surgery;
- have had any recent illness or infections;
- have any other medical conditions or complaints.

Before using the massager, consult your doctor regarding any unexplained pain, if you are undergoing medical treatment or using medical equipment. Discontinue using the Proberg MG-999 if you experience pain, discomfort, or any abnormal symptoms.

What's in the box

1	PROBERG MG-999 Percussion massager
7	Silicone massage attachments
1	Power adapter
1	User manual
1	Case

Parts of the product

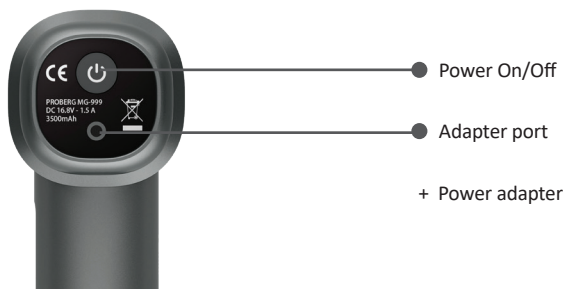Spherical attachment

Available in 3 sizes: large (diameter 55 mm), medium (diameter 50 mm), small (diameter 43 mm). Suitable for arms, waist, back, thighs, calves and other large muscle groups.

Flat attachment

Suitable for various muscle groups.

U-shaped attachment

Used for targeted massage of the deep muscles along the spine, neck and Achilles tendon.

Cylindrical attachment

Suitable for deep tissue and trigger points such as neck, palms, soles.

Dual-ball attachment:

Ideal for targeting larger muscle groups and providing a soothing massage experience. Can help relieve tension in the back, thighs and calves.

DO NOT use any of the attachments on your head!

Charging the Proberg MG-999

1. When using for the first time, please fully charge the battery for approximately 3.5 hours. The battery run time when fully charged is 3 to 8 hours depending on the speed level you choose.
2. To charge the device turn it off, connect the power adapter, the power bar indicator light is green and the percentage indicator light is green and flashes in sequence.
3. Once the device is fully charged, the power bar indicator light is solid green and the light is no longer flashing.
4. The device cannot be used while charging!

Instructions for use for the Proberg MG-999

1. When using for the first time, please fully charge the battery for approximately 3.5 hours.
2. Choose a suitable massage attachment and screw it into the socket on the front of the device.
3. Press and hold the power button for 2 seconds to turn on the device. The product enters the lowest speed mode. The device features 8 speed levels. Use the speed selector on the top of the device to adjust the speed by slightly touching the "+" or "-" button. The speed indicator on the back of the device will illuminate the speed chosen. If you have not used this type of product before it is recommended that you start on level 1 speed and gradually increase its intensity. Once you have finished with the massage, press and hold the power button for 2 seconds to turn off the device.
4. The device will automatically switch off after 15 minutes of continuous use.
5. The battery level indicator on the back of the device shows the current battery charge level. The power bar corresponds to the percentage ratio (25% 50% 75% 100%) indicating different power range. The bars on the battery will show accordingly the battery percentage.
6. When the battery is low, the battery light flashes 5 times and the device automatically turns off.

! Overload protection: If too much pressure is applied to the attachment, the device will shut off automatically. Release the attachment and the device will resume work.

Technical specifications for the Proberg MG-999

Battery capacity	3500 mAh/14.4V
Battery type	SAMSUNG Rechargeable
Charging port	DC - 16.8 V / 1.5A
Charging time	Approx. 3.5 hours
Battery run time	Approx. 3-8 hours
Input Voltage	DC 16.8V
Rated power	33W
Speed	1500-2900 percussions per minute
Amplitude	12 mm
Net Weight	0.98 kg (excluding massage heads)
Gross Weight	2.5 kg
Product size	18.7 x 11.4 x 5.8 cm
Massage attachments	7
Colour	Grey
Automatic switch off	After approx. 15 minutes

Cleaning and maintenance

- Make sure the device is turned off, the cable unplugged and allow it to cool before cleaning. DO NOT use abrasive cleaners, such as benzene or thinner to clean the product, as this may cause the colour of the product to fade.
- DO NOT immerse the device in water or any other liquids during the cleaning or at any other times.
- During the cleaning, prevent water or other liquids from entering the device to avoid short circuit and electric shock hazards during use.
- Remove any small marks with a damp cloth.
- Store the device in its original case in a cool, dry place.
- Do not place any other heavy weight/objects over the device to avoid malfunction and/or damage.

Troubleshooting

PROBLEM	POSSIBLE CAUSE	SOLUTION
The device does not turn on.	Power is off.	Press and hold the power button for 2 seconds to turn on the device.
	Low battery	Fully charge the device and then turn it on.
The device automatically switches off.	Overheat protection: If the device has overheated, it will shut off automatically.	Leave the device for half an hour to cool down and then try turning it on again.
	The device automatically shuts off after 15 minutes of continuous use. (Timer)	Press and hold the power button for 2 seconds to turn on the device.
	Overload protection: If too much pressure is applied to the attachment, the device will shut off automatically.	Release the attachment and the device will resume work.
If you encounter any other problems, please contact our customer support at de.info@donnerberg.net		

Warranty policy

The company Donnerberg S.G., Von-Werthern-Str. 4a, 96487 Dörfles-Esbach, Germany grants a valid manufacturer's warranty for the territory of the European Union and the civil-law defect liability. The warranty starts from the time the device was delivered to the customer. Battery is excluded from this warranty. It is not extended by a previous warranty case. The warranty is granted only to the buyer and is not transferable to third parties.

In addition to the statutory liability for defects, we also repair or replace equipment whose defects are not initial defects unless they are caused by improper use or are insignificant. The choice between replacement and repair is ours. The warranty is only granted if the defect is notified immediately upon its occurrence and is reported. You must include the claim form and the copy of your purchase receipt as the proof of warranty. Please note that we will not process warranty claims without the claim form. Upon notifying our customer service team, we will send you a claim form.

Electric and electronic devices

Information for private households

The Waste Electrical and Electronic Equipment Regulations 2018 (WEEE) in the United Kingdom contains several requirements concerning the handling of electric and electronic devices as waste. In the following paragraphs, we list the most important requirements for you.

1. Separated recycling of „old devices“

Electric and electronic devices that have become waste are called old devices. Owners of such old devices have to bring these to recycling facilities, where these are handled and recorded separately from urban solid waste. Old devices must not end up in the household refuse. In fact they have to be collected in special collecting and recycling systems.

2. Possibility for disposal of „old devices“

Owners of old devices from private households can dispose these at the public waste disposal authorities or bring them to producer or seller related disposal locations. You can find these disposal locations online: <http://www.recycle-more.co.uk>

3. Privacy information

Old devices often contain private data, especially appliances from data processing and communication equipment, as smartphones or hard discs. In your own interest, please arrange to eliminate all personal and private data from your disposed equipment.

4. Meaning of the symbol „crossed out waste container“

The labelled symbol showing the „crossed out waste container“ on electric and electronic equipment indicates that the concerning equipment has to be disposed separately from urban solid waste.

5. The Waste Electrical and Electronic Equipment Regulations 2018 (WEEE)

Following the implementation of the Directive 2012/19/EU and its British transportation WEEE Regulations 2018, we as a producer of EEE are registered in the Public register of Electrical and Electronic Equipment Producers. The register is held in England and Wales by the Environment Agency.

Guidance for safe removal of batteries & accumulators

This electrical device contains the following battery:

Battery type	Chemical composition
ES18650-4S1P 3500mAh 14.4V	Lithium Ion Battery

Please do not remove the battery yourself within the warranty period. This invalidates the manufacturer's warranty. Please contact our customer service at de.info@donnerberg.net.

At the end of the product's service life, you can dispose of the battery yourself for recycling as described:

Warning: Make sure that the battery is completely discharged and that the device is switched off and the charging cable unplugged.

1. Utilize tweezers to carefully peel off the bottom sticker surrounding the round head charging port.
2. With a Phillips screwdriver, undo the two screws securing the bottom cover.
3. Disconnect all terminal wires from the PCBA terminal block, totaling three terminals.
4. Employ tweezers to extract the silicone pads located on both sides of the handle.
5. Utilize a Phillips screwdriver to remove the four screws fastening the handle in place.
6. Detach the handle from the main body of the device.
7. Using a Phillips screwdriver, remove the four screws situated on both sides of the handle.
8. Gently push the battery out from within the handle. The battery and the device can now be disposed of separately.

Battery safety information

The following precautions apply to all products that use rechargeable batteries, including lithium ion battery. Improper battery use may result in serious injury, death, property damage, or damage to the product as a result of battery fluid leakage, fire, overheating or explosion. Released battery fluid is corrosive and may be toxic. It can cause skin and eye burns, and is harmful if swallowed. If a battery leaks, remove all batteries by reversing the installation steps provided for this battery, being careful to keep the leaked fluid from touching your skin or clothes. If fluid battery contacts skin or clothes, flush skin with water immediately.

To reduce risk of injury and proper performance of your battery:

- Keep batteries out of the reach of children.
- If a battery has been swallowed seek medical attention immediately.
- Do not dismantle, open or shred battery pack or cells.
- Do not subject batteries to mechanical shock.
- Only replace with batteries of the same type, never use different types of batteries together or used batteries with new ones. Do not mix different size, type, manufacture, or capacity within the massager.
- Remove the batteries if the massager will not be used for an extended period of time.

- Do not short-circuit the battery. Do not store the battery in a box or a drawer where it may short-circuit each other or be short-circuited by other metal objects.
- Do not overcharge batteries! There is a risk of explosion!
- Do not short-circuit battery, as there is danger of explosion.
- Do not throw the battery into fire, as there is danger of explosion.
- Do not throw used batteries in the household waste, but in the hazardous waste or take them to designated battery collection points.

Customer service

We have a reliable customer service team that is committed to providing solutions to any issue you might have. If you have any questions or concerns about the products, warranty or exchange do not hesitate to contact us at:

E-mail	de.info@donnerberg.net
Telephone	+ 49 89 215 400 980 / + 49 95 619 769 199
Donnerberg S.G. Von-Werthern-Str. 4a 96487 Dörfles-Esbach Germany	

Translations of the User Manual into other languages can be found at www.donnerberg.net.

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