PREMIUM AKKU

NM-090

SAMSUNG
2600mAh
AUFLADBARE
BATTERIE

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ENGLISH

Contents

Introduction	15
Traditional shiatsu massage.....	15
Intended use	15
Product description	15
Parts of the product.....	16
Samsung rechargeable battery.....	16
Product functions.....	17
Product positions	18
Safety precautions	19
Specific group of people	20
Contents included	20
Technical specifications.....	21
Maintenance and storage.....	21
Electric and electronic devices	22
Guidance for safe disposal of batteries & accumulators.....	23
Battery safety information.....	23
Certificates and recommendations	24
Customer service	25
Warranty policy.....	25

Introduction

Dear customer,

Thank you for choosing Donnerberg Original Premium neck and shoulder massager with rechargeable battery.

Please read the manual instructions carefully and save it for later use, be sure to make them accessible to the other users. Please follow the instruction and use the product as described.

Donnerberg Team

Traditional shiatsu massage

Donnerberg Original Premium massager with rechargeable battery is based on the principles of the Shiatsu massage that originated in Japan. This massage supports the body's natural ability to self-heal. The word Shiatsu means 'finger pressure' in Japanese. Shiatsu uses pressure, touch and manipulative techniques to bring the body and its inner energy paths into balance. Shiatsu massage can be used to reduce stress and promote overall health and well-being.

Intended use




The unit is only intended for domestic /private use, not for commercial use. It is designed for massaging the following body parts: neck, back, shoulders, belly, legs and arms. If you are unsure whether the neck and shoulder massager is suitable for you, please consult your doctor.

Product description

Sit in a comfortable chair, carefully place the massager over your shoulders and switch it on. Put the arms through the loops to facilitate the positioning of the massager and to control the intensity of the massage. Place the massager on the desired areas or parts of the body and adjust the intensity of the massage by pulling harder or softer on the loops.

Each massage cycle lasts 15 minutes but can be shortened on demand. When the 15 minute cycle has elapsed, the massager will automatically switch off. This device must be turned off when not in use.

Parts of the product

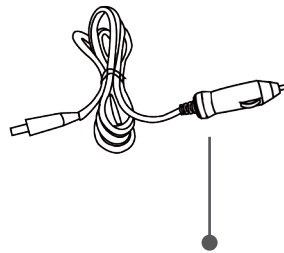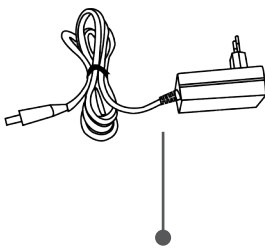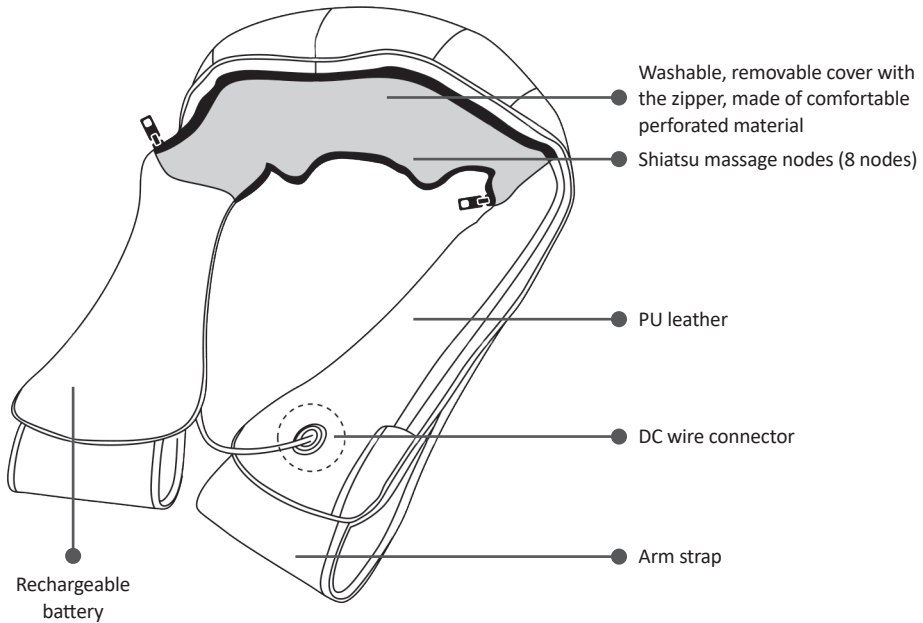Samsung rechargeable battery

Rechargeable massager with SAMSUNG battery:

This massager has a rechargeable battery, which allows its autonomous use without plugging into the outlet, if the battery is fully charged. You can use this device independently with the battery for approximately 80 minutes if the battery is completely full. Complete charging of the battery lasts around 3 hours. The battery charges when the massager is plugged into a socket.

A solid blue light indicates that the battery is charging. Once the battery is fully charged, the blue light will turn off. Samsung rechargeable battery is located on the arm strap opposite of the console.

Product functions

1. Power on/off

- Press this button to turn on the device. When the battery is the main power supply, hold this button for three seconds in order to switch on the massager.
- This device is equipped with an infrared heating function which automatically turns on by pressing the “Power” button along with the massage function and rotation.
- Press the “Infrared heat” button to turn off the heating function.

2. Direction

- Press the “Direction” button to change the rotating direction.

3. Speed

- Choose between 3 massage speed modes. The lowest speed is preset by default when the massage device is switched on.
- Pressing the “Speed” button increases the speed of the massage.
 - Press 1x for the medium speed.
 - Press 2x for the highest speed.

4. Infrared heat

- You can choose between two levels of infrared heat:
 - When the device is switched on, the 2nd level of the Infrared heat is on.
 - Press 1x to switch to the 1st level of infrared heat.
 - Press 2x to turn off the infrared heat.
- You can also use infrared heat independently of the massage.
- To activate only infrared heat, switch off the device and press the “Infrared heat” button.






5. Vibration

- Turn on the additional vibration function by pressing the “Vibration” button in case you would like more intensive massage.
 - Press 1x to turn on vibration. Default vibration mode - continuous vibration.
 - Press 2x to change the speed and mode of vibration. Vibration lasts for 5 seconds, then pauses for 2 seconds repetitively.
-
- The vibration function is only recommended on the shoulders, back and leg area.
 - The massage should feel pleasant and relaxing at all times. If the massage is painful or at any time uncomfortable change the pressure using the straps.

Product positions

Please apply on body parts shown in the pictures below.

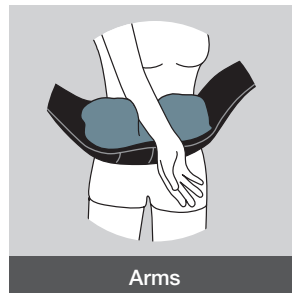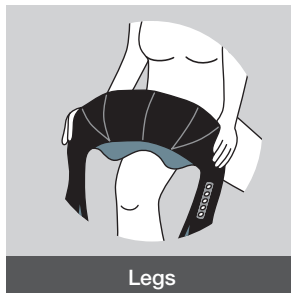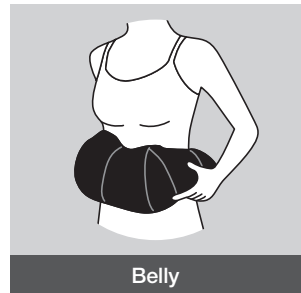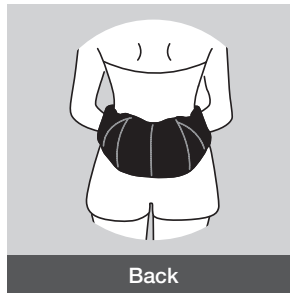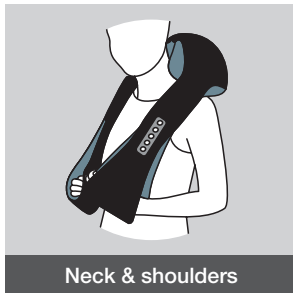Safety precautions

To prevent the risk of fire, electric shock, or other injury, please follow the instructions:

1. Before connecting the device to your power supply, please ensure that the supply voltage stated on the rating plates is compatible with your mains supply. Only insert the power supply unit into a socket when the appliance is switched off.
2. Do not cover or use the unit under a blanket when operating. Do not use while bathing or showering. The unit must be kept dry and stored in a dry and safe place. The power socket, cable, and appliance must be at appropriate distance from the heat, hot surfaces, moistures and liquids. Connect the unit only in that way, so that the power plug is freely accessible. The power plug should be freely accessible when inserted into the socket. The cable must not be kinked, trapped or twisted.
3. Change the covers only when the device is turned off, the power supply unit unplugged and the battery fully discharged.
4. To use this product safely and correctly, read this manual thoroughly and operate the device as instructed. Never operate this appliance if it has a damaged cord or plug or if it is not working properly. In case the device is not working or has a fault, send it to the service centre for examination. Do not plug in the device in the damaged sockets. If the device has been dropped into water, do not pick it up or touch it. Pull out the power plug from the socket immediately.

Do not operate where aerosol (spray) products are being used or where oxygen is being administered. After you have finished using the unit turn off the power. Before unplugging, the device must be turned off. Do not obstruct with objects or thrust.

Please ensure that the polythene packing is kept away from the reach of children!
Risk of suffocation and choking hazard!

Do not use the massager while driving!

Children are not allowed to play with the device!

Specific group of people

- This appliance can be used by children aged 8 and over and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved.
 - Children are not allowed to play with the appliance. Cleaning and user maintenance must not be done by children without adult supervision.
 - The appliance has a heated surface. Persons sensitive to heat must be careful when using the appliance.
 - Do not use this device to supplement or replace medical treatment.
 - Chronic complaints and symptoms may be made worse.
 - Do not use the appliance near your eyes or near any other sensitive areas.
 - Take care when using the heating function. Do not use if you have circulatory problems or on areas sensitive to heat.
 - If you experience any pain or discomfort during the massage, stop immediately and consult your doctor.
 - Before using the massager, consult your doctor regarding any unexplained pain, if you are undergoing medical treatment or using medical equipment.
- The appliance is only to be used with the power supply unit provided with the appliance.
 - This appliance contains a non user-replaceable battery.
 - This appliance contains a battery that is only replaceable by skilled persons or authorized service centres.
 - For the purposes of recharging the battery, only use the detachable power supply unit provided with the appliance.
- The model number of the detachable supply unit is IVP1280-2000G.

You should refrain from using the Shiatsu massager, or should consult your doctor before using if:

- You are pregnant,
- You have a cardiac pacemaker, artificial joints or electronic implants,
- You suffer from one or more of the following illnesses and complaints: Circulation disorders, varicose veins, open wounds, bruises, broken skin, inflammation of the veins.

Contents included


1	Donnerberg Original Premium Neck & Shoulder massager with rechargeable battery NM-090
2	Removable covers in different colours
1	Power adapter (length approx 2.8m)
1	Car charger
1	Soft cloth cover for sensitive skin
1	Protective custody - environmentally friendly
1	User manual
1	Samsung 2600mAh Battery rechargeable

Technical specifications

Product name	PREMIUM AKKU NM-090
Product model	NM-090 Black
Voltage supply	Output DC 12.8V 2A 25.6W
Input voltage	100-240 V, 50/60 Hz
Automatic switch off	After approx. 15 minutes
Net Weight	Approx. 1,98 kg
Gross Weight	Approx. 2,65 kg
Product dimensions	Approx. 50 x 43 x 17,6 cm

Maintenance and storage

- Before cleaning, make sure the switch is turned off and power plug is disconnected from the outlet. Do not use abrasive cleaners such as benzene or thinner to clean the product, otherwise the colour of the product may fade. Do not immerse the appliance in water or the other liquids during the cleaning! Small marks can be removed by using a damp cloth. Please remember to clean only the top side of the device.
- You may only clean the massager yourself. If the device malfunctions please do not repair it yourself as this leads to automatic expiry of all guarantee rights.
- Thanks to the integrated zipper, the cover can be easily removed and washed in the washing machine at 30° C. Simple and easy, just unzip the zipper and replace the cover with another one.

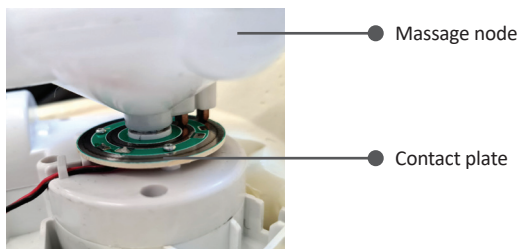- We would kindly advise you to clean your massager regularly by taking off the cover and removing dirt (hair, dust) that gathered between the massage node and the contact plate with cotton buds. Regular removal of contamination can prevent malfunctions of certain components of infrared heat and massage.

Electric and electronic devices

Information for private households

The Waste Electrical and Electronic Equipment Regulations 2018 (WEEE) in the United Kingdom contains several requirements concerning the handling of electric and electronic devices as waste. In the following paragraphs, we list the most important requirements for you.

1. Separated recycling of „old devices“

Electric and electronic devices that have become waste are called old devices. Owners of such old devices have to bring these to recycling facilities, where these are handled and recorded separately from urban solid waste. Old devices must not end up in the household refuse. In fact they have to be collected in special collecting and recycling systems.

2. Possibility for disposal of „old devices“

Owners of old devices from private households can dispose these at the public waste disposal authorities or bring them to producer or seller related disposal locations. You can find these disposal locations online: <http://www.recycle-more.co.uk>

3. Privacy information

Old devices often contain private data, especially appliances from data processing and communication equipment, as smartphones or hard discs. In your own interest, please arrange to eliminate all personal and private data from your disposed equipment.

4. Meaning of the symbol „crossed out waste container“

The labelled symbol showing the „crossed out waste container“ on electric and electronic equipment indicates that the concerning equipment has to be disposed separately from urban solid waste.

5. The Waste Electrical and Electronic Equipment Regulations 2018 (WEEE)

Following the implementation of the Directive 2012/19/EU and its British transportation WEEE Regulations 2018, we as a producer of EEE are registered in the Public register of Electrical and Electronic Equipment Producers. The register is held in England and Wales by the Environment Agency.

Guidance for safe disposal of batteries & accumulators

This electrical device contains the following battery:

Battery type	Chemical composition
FBC 18650 2600mAh 11.1V	Lithium Ion Battery

Please do not remove the battery yourself within the warranty period.

This invalidates the manufacturer's warranty. Please contact our customer service at de.info@donnerberg.net.

At the end of the product's service life, you can dispose of the battery yourself for recycling as described:

Warning: Make sure that the battery is completely discharged and that the device is switched off and the charging cable unplugged.

- Open the outer cover with the zipper.
- You will find the battery box in the left arm loop.
- Remove the strips around the battery and disconnect the wire connector.
- Carefully remove the battery or accumulator.
- The battery or accumulator and the device can now be disposed of separately.

Battery safety information

The following precautions apply to all products that use rechargeable batteries, including lithium ion battery. Improper battery use may result in serious injury, death, property damage, or damage to the product as a result of battery fluid leakage, fire, overheating or explosion. Released battery fluid is corrosive and may be toxic. It can cause skin and eye burns, and is harmful if swallowed. If a battery leaks, remove all batteries by reversing the installation steps provided for this battery, being careful to keep the leaked fluid from touching your skin or clothes. If fluid battery contacts skin or clothes, flush skin with water immediately. Before inserting new batteries, thoroughly clean the battery compartment with a dry cloth.

To reduce risk of injury and proper performance of your battery:

- Keep batteries out of the reach of children.
- If a battery has been swallowed seek medical attention immediately.
- Do not dismantle, open or shred battery pack or cells.
- Do not subject batteries to mechanical shock.
- Only replace with batteries of the same type, never use different types of batteries together or used batteries with new ones. Do not mix different size, type, manufacture, or capacity within the massager.
- Observe the plus (+) and the minus (-) marks on the battery to ensure the correct use
- Remove the batteries if the massager will not be used for an extended period of time
- Do not short-circuit the battery. Do not store the battery in a box or a drawer where it may short-circuit each other or be short-circuited by other metal objects.
- Do not overcharge batteries! There is a risk of explosion!
- Do not short-circuit battery, as there is danger of explosion.
- Do not throw the battery into fire, as there is danger of explosion.
- Do not throw used batteries in the household waste, but in the hazardous waste or take them to designated battery collection points.

Charging protection scheme:

1. Overcharge protection: The battery can be charged up to 12.6V. The overcharge protection is protected by software. The output voltage of the external switching power supply is controlled at 12.6V~12.9V, which also controls the maximum voltage of the battery, which is 12.6V.
2. Discharge protection: When the product is powered by battery, the battery discharge will be stopped as low as 8.6V. After the battery is discharged, you need to re-insert the power supply to work. This protection is set by software.
3. Charging time protection: After 3h, charging will immediately stop
4. Charging current protection: charging current is controlled within 2A.
5. Short circuit protection: When the battery is connected between positive and negative, the battery built-in protection board will cut off the power.
6. Battery internal protection: The battery has PCB protection; overcharge setting is set to 12.6V and the over discharge setting is set to 8.25V.

Certificates and recommendations

 	<p>Recommended by the German Olympic Sports Confederation</p>
	<p>TÜV SÜD (Technical Inspection Association)-the product was tested for electrical safety, hazardous material components, and unwanted noise production. Certificate No.: Z2 091698 0007 Rev.02</p>
	<p>RoHS (Restriction of Hazardous Substances) - restricts the use of specific hazardous materials found in electrical and electronic products 2011/65/EU</p>
	<p>CE (Conformite Europeenne)- indicates that the product conforms with the relevant EU directives regarding health and safety or environmental protection (765/2008)</p>

Customer service

We have a reliable customer service team that is committed to providing solutions to any issue you might have. If you have any questions or concerns about the products, warranty or exchange do not hesitate to contact us at:

E-mail	de.info@donnerberg.net
Telephone	+ 49 89 215 400 980 / + 49 95 619 769 199
Company Donnerberg Von-Werthern-Str. 4a 96487 Dörfles-Esbach Germany	

Warranty policy

The company Donnerberg, Mr. Svetozar Grbic, Von-Werthern-Str. 4a, 96487 Dörfles-Esbach, Germany grants a valid manufacturer's warranty for the territory of the European Union and the civil-law defect liability. The warranty starts from the time the device was delivered to the customer. Battery is excluded from this warranty. It is not extended by a previous warranty case. The warranty is granted only to the buyer and is not transferable to third parties.

In addition to the statutory liability for defects, we also repair or replace equipment whose defects are not initial defects unless they are caused by improper use or are insignificant. The choice between replacement and repair is ours. The warranty is only granted if the defect is notified immediately upon its occurrence and is reported. You must include the claim form and the copy of your purchase receipt as the proof of warranty. Please note that we will not process warranty claims without this claim form.

You may report your warranty claim to one of the following:

**Company Donnerberg
Von-Werthern-Str. 4a
96487 Dörfles-Esbach
Germany**

e-mail: de.info@donnerberg.net

or telephone: + 49 89 215 400 980 / + 49 95 619 769 199

Upon your notification, our customer service team will send you a claim form.

You must include the claim form and the copy of your purchase receipt as the proof of warranty. Please note that we will not process warranty claims without this claim form.

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